

Welch Allyn PC-Based SpiroPerfect Spirometer

Technical Specifications

SpiroPerfect Tests	Computer based full diagnostic spirometer FVC, SVC, MVV, Pre-Post BD	Quality Checks	Effort Acceptability and Test Reproducibility checks Audio and visual incentive for assistance in coaching patients
Sensor Type	Pneumotach	Connectivity	Part of the CardioPerfect Diagnostic Suite Export compatible with most electronic medical records programs Available in multiuser network Telemedicine option for e-mail transfer
Power Equipment	None, obtained from USB port		
Accuracy	Meets or exceeds 1994 ATS standard		
Reproducibility	Meets or exceeds 1994 ATS standard		
Volume Range	0-14L		
Flow Range	±14L/sec		
Software	Windows 2000 and XP compatible		
Predictive Norms	Berglund 1963, Crapo 1981, Composite, Dockery 1983, ECCS / Quanjer 1993, Falaschetti 1995, Gulsvik 2001, Hedenström 1986, Hsu 1979, Knudson 1976, Koillinen 1998, Kory 1961, Morris 1971, NHANES III 1999 [Default], Polgar 1971, Schoenberg 1978, Solymar 1980, Viljanen 1981, Zapletal 1969 <small>Additional predictive norms can be added upon customer request</small>	Minimum Computer Requirements (Recommended)	
Interpretation	1991 ATS Interpretation Standards. Automatic can be disabled. Manual Interpretation available. Lung Age calculation	Processor	Intel® Pentium® III (Pentium IV or higher)
Reports	Multiple report templates including FVC, SVC, MVV, measurement, and trend	Hard Drive	2GB (4GB or more)
Incentive Graphic Parameters	Fireman FVC, FIVC, FIV1, FIV1%, FEV0.5, FEV1, FEV2, FEV3, FEV5, FEV6, FEV1/FEV6, FEV0.5%, FEV1%, FEV2%, FEV3%, FEV5%, FEV6%, PEF, FEF25, FEF50, FEF75, FEF0.2-1.2, FEF25-75, FEF75-85, PIF, FIF50, FEF50/FIF50, FET SVC, ERV, IRV, VT, IC, BF, MV, Tin, Tex, Tin/Tex	RAM Memory	256MB (512MB or higher)
		Monitor Resolution	(1024 x 768 true color)
		Operating Systems	Microsoft® Windows® 2000 or XP
		CD-ROM	Required for installation
		Printers	Color printer recommended. Compatible with inkjet and laser printers. Laser printers require a minimum of 2MB of memory
		Performance Standards	SpiroPerfect meets or exceeds the following standards: ATS—American Thoracic Society, Standardization of Spirometry, 1994 Update ERS—European Respiratory Society NIOSH—National Institute for Occupational Safety and Health OSHA—Occupational Safety and Health Administration SSD—Social Security Disability ECCS—European Community for Coal & Steel

Ordering Information

SPIRO-S	CPWS PC Spirometer Module, with 3L calibration syringe. Includes: sensor, pressure tubing, 4 flow transducers, electronic operator's manual, latest CPWS SW CD, 2 training posters, and quick reference card
SPIRO-NS	CPWS PC Spirometer Module, without 3L calibration syringe. Includes sensor, pressure tubing, 4 flow transducers, electronic operator's manual, latest CPWS SW CD, 2 training posters, and quick reference card
703418	Disposable Flow Transducer, 25 pk, CPWS, CP 200
703419	Disposable Flow Transducer, 100 pk, CPWS, CP 200
703415	Pressure Tubing, CPWS, CP 200, 2m
58550-0000	Nose Clip (1 each)
703480	Calibration Syringe (3 liters), CPWS, CP 200

To learn more about the Welch Allyn PC-Based SpiroPerfect Spirometer or to request a demonstration, contact your Welch Allyn representative. You'll find a solution that not only improves patient care, but provides priceless peace of mind.


4341 State Street Road, PO Box 220, Skaneateles Falls, NY 13153-0220 USA
(p)800.535.6663 (f)315.685.3361 www.welchallyn.com

© 2005 Welch Allyn SM2846 Rev. B

WelchAllyn

Advancing Frontline Care™

Cardiopulmonary Solutions

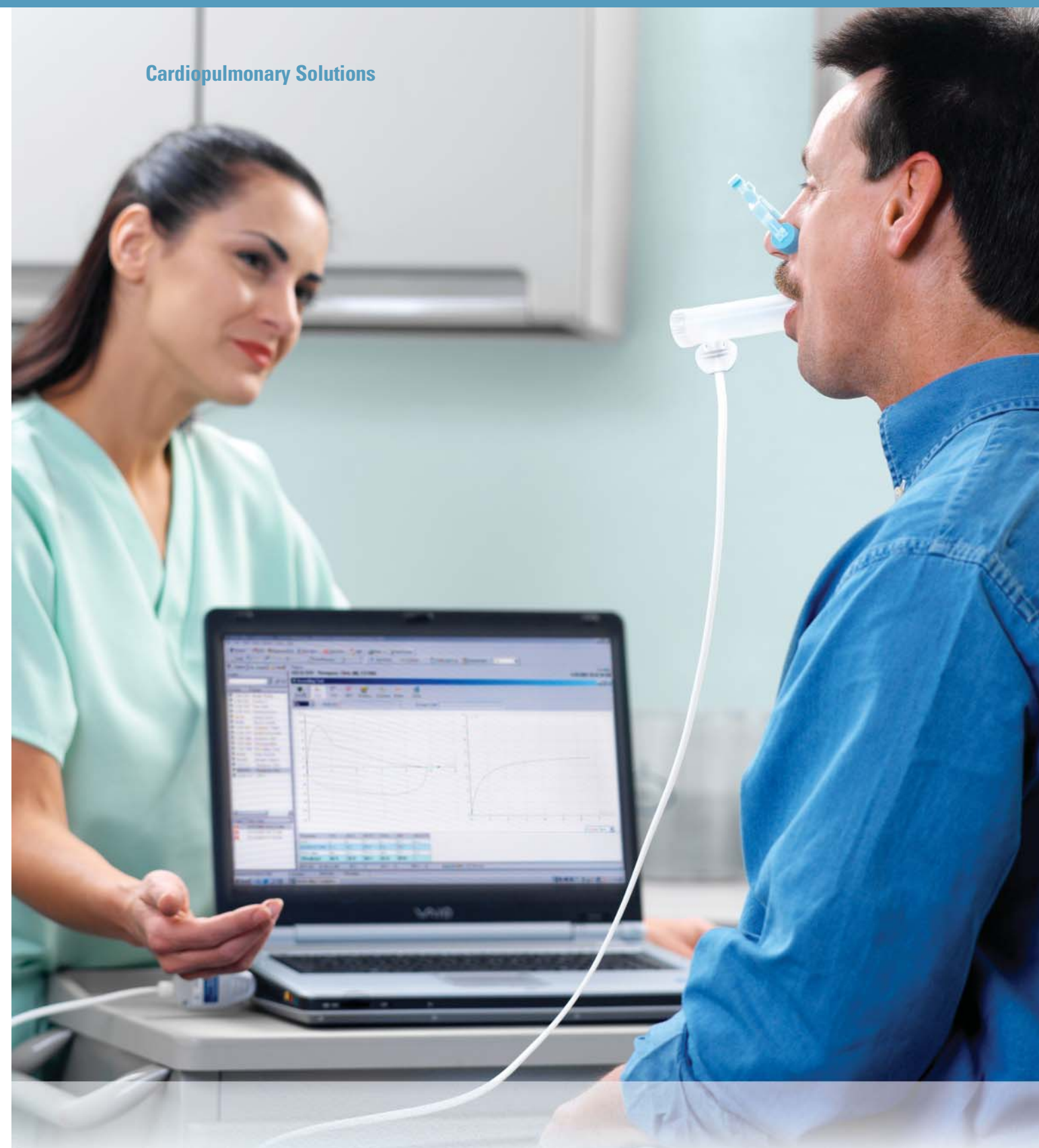

Welch Allyn PC-Based SpiroPerfect™ Spirometer
Precise information. Easier diagnoses. Better patient care.


WelchAllyn

Advancing Frontline Care™

Accurate, integrated patient information at your fingertips

In seconds, you and your staff can see and compare accurate results from pulmonary function tests with the friendly, feather-weight 142g SpiroPerfect. You can clearly evaluate pulmonary obstruction or restriction, increase success rates of patients who are undergoing smoking cessation, show progress, and help pediatric patients comply with test procedures.

- Fast, reliable Pre-Post bronchodilator FVC, SVC and MVV testing with automatic interpretation and reversibility analysis
 - > View real-time Flow/Volume Loop and Volume/Time graphs
 - > Instantly verify quality and variability of proper test performance
 - > Motivate children with incentive graphic—for pediatric patient coaching
 - > Reduce risk of cross contamination with single-use disposable transducers
 - > Select single and multiple flow calibration protocols
- Easily monitor, review, store and print on your existing computer with the Welch Allyn CardioPerfect™ Workstation
 - > Choose from multiple predicted norms and report formats
 - > Compare previous tests with the Trending tool
 - > Make available to multiple staff via computer network
 - > Interface with Electronic Medical Records or patient management system for onsite or offsite review of key data. Telemedicine option is available for e-mail transfer.

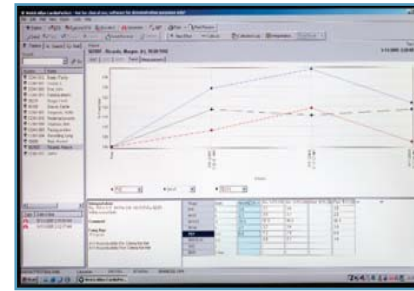

Compare tests done on different days and clearly see progress. View three parameters at once as a percentage of predictive values or of absolute value.




Enjoy the convenience of our accurate patented **Disposable Flow Transducers**. Designed for use by a single adult or child, the breakaway feature minimizes risk of cross contamination. No sterilization, no cleanup—just safety and savings.



Coach, motivate and make testing more fun. Because children understand how to help the fireman blow out the fire, the **Incentive Screen** helps kids understand your goals and comply.

Get Connected to Your EMR through Welch Allyn's CardioPerfect™ Workstation Software

Integrated. Customizable. Intuitive.

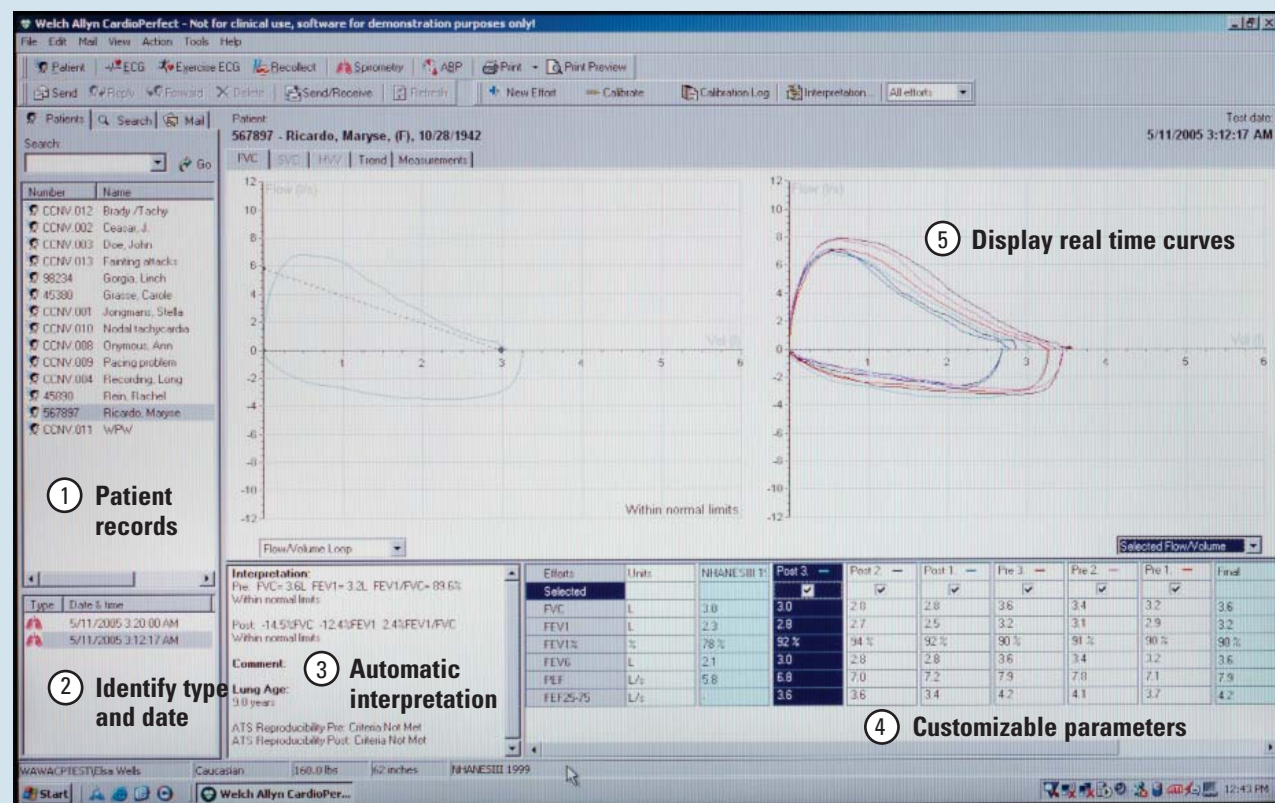

- 1 Store, retrieve, print or edit **patient records** using the Welch Allyn CardioPerfect Workstation database.
- 2 Easily **identify** all diagnostic test types and dates each procedure was performed.
- 3 **Automatic interpretation** follows ATS Interpretative Strategies and test reproducibility criteria.
- 4 See a quick snapshot of patient results. Select from 6 parameters to show real-time testing and summary. View measurement tables with 19 fully **customizable parameters**.
- 5 **Display real time curves** as Flow/Volume, Volume/Time, Flow/Time or Pre, Post. All effort or best effort can be displayed on the same graph for easy comparison.



Welch Allyn PC-Based SpiroPerfect™ Spirometer



Welch Allyn CardioPerfect™ Workstation Software



Connects to Electronic Medical Records Systems

The Welch Allyn PC-Based SpiroPerfect™ Spirometer downloads information to the Welch Allyn CardioPerfect™ Workstation Software, which provides paperless efficiency as well as connectivity to many leading EMR systems. This software puts patient records from a variety of Welch Allyn diagnostic instruments at your fingertips, giving you the ability to store, access, review, compare and edit diagnostic test data from any computer on your network.

- Retrieve records with the touch of a key, and review and edit right on the screen
- Store thousands of records with ease
- Access information from your exam rooms, your office, or at home
- Fax or e-mail reports directly for consults
- Add additional cardiopulmonary diagnostic testing devices like ECG, ABPM, or Holter as your practice grows
- Copy and paste spiromographs into Microsoft® Word
- Integrate with other patient information systems like billing
- Uses passwords to assign/limit tasks so that data integrity is maintained